

Product Catalog Lab Range

Shiva Products

ISO 13485 : 2016 (TUV SUD Certified)

Proudly Exporting to 40+ Countries Worldwide

2021



INDEX



Cold Cure Powder / Liquid	5
Cold Mould Seal	5
Heat Cure Powder / Liquid	6
R2 Alloy for PFM	
Ni-Cr Alloy for C/B	
Die Hardner	8
Die Spacer Thinner Set	8
Die Separator	
Reinforcement Mesh	
Modelling Wax	10
Shellac Base Plates	10
Dental Stone Plasters	11
Ortho White Stone Plaster	11
Model Stone Plaster	12
Milling Wax Blanks	12
Sprue / Wire wax	13
Lingual Bar Wax	13
Utility Wax	14
Grids Wax	14
Retention Waxes	15
Clasp Waxes	15
Inlay / Onlay Pattern Wax	16
Die Stone Plaster (Pink)	17
Casting Wax Sheets	17

INSTITUTIONAL RANGE

Carving Wax / Wax Blocks	18
Beading Wax	18
Green Tracing Sticks	19
In l ay Sticks	19
Sticky Sticks	19
Spacer Wax	20
Boxing Wax	20



COMPANY PROFILE

An ISO 13485 : 2016 accredited company catering a wide range of products such as dental material, dental products, dental alloy, jewellery wax, etc...

An Overview

Shiva Products, established in the year 1998,an eminent company is a manufacturer for Dental Products and other specialty waxes. With more than two decades of experience in this industry, we are proud to be recognized as a renowned manufacturer and exporter of premium quality Dental Goods, Jewellery Waxes, Investment Casting Waxes, Modelling Waxes, Dental Alloys, Dental Instruments, etc., which are specifically made according to international quality parameters. We are a team of experienced professionals working with a common goal of quality production, thereby, winning the trust and faith of thousands of national and international clients.

The Company has ventured into the field of dentistry since 1998 with 2 products in hand. With Team work & focus on RRR strategy (Reliable Responsible & Reasonable), we currently have more than 80 Products catalogued and many in pipe line. Company has become on of the leading manufacturer in Domestic & well-known in International Market. The Strong relationship and the trust of consumers, trade and the dental profession built over two decade of operations in India has made 'MAARC' a trusted Brand.

'MAARC' products have the reputation of "RIGHT QUALITY RIGHT PRICE" in India as well as International market. Our ISO 13485:2016 manufacturing standards helps us to deliver quality dentistry products at affordable rates to masses. Our goal is to build confident relationship with our customers.

We also cater to manufacturers in neutral packing and contract manufacturing can also be done.

DENTAL SEGMENT:

Clientale:

01. Domestic: We are currently selling to approx. 1000-1200 Distributor / Dealers (many states of India), approx. 500-600 Laboratories within Mumbai, 4000-5000 dentists are been catered directly / through distributor/dealer network through our company within Mumbai + Institutional supply + government Tender supply.

02. Exports - Currently supplying to approx. 40 countries across various continents

Certifications: Shiva Products has several certifications as under.

- a) EN ISO 13485: 2016
- b) Class I device products are all CE certified
- c) EAR Obelis SA -Belgium
- d) Registration with various Chamber of commerce in Mumbai for smooth functioning of export orders.



COMPANY PROFILE

Vision

We are led by global vision and drive to manufacture various dental material comparable to international standards so as to offer optimal customer satisfaction.

With a goal to establish strong relationships with more and more customers, we offer our Products at extremely affordable price. It is also a continuous effort on our part to deliver the orders within the stipulated time.

Assured Quality

We take utmost care of quality while producing our goods, which meet international quality norms in all respects. Shiva Products has been certified with an ISO 13485: 2016 certification, which is a testimony in itself to our quality production. Each operation is executed under the strict supervision of trained professionals in accordance with standardized production norms and methods, thus leaving no stone unturned in ensuring premium quality.

Infrastructure

Shiva Products has a state-of-the-art manufacturing unit, where latest manufacturing machines and equipment contribute substantially in manufacturing global quality dental products. Our manpower includes engineers, dentists, technicians, designers, R&D team, QC, marketing professionals, quality personnel in addition to some semi -skilled workers. We work as a team to achieve our vision.

In-process Quality Analysis: Our manufacturing unit has its own lab and checking facilities wherein quality analysis is been done on daily basis for RM, in-process and final product before its sent to the packing zone.

The most advantageous and beneficial technique of production supplied by our R&D team give us an edge over our competitors. Our warehouse is a boon, as it facilitates safe storage of bulk products to meet all types of market demands.

Future Plans

Our Company want to take advantage of the booming export market to add to its market share and not lose the opportunity to get an insight into the behavior of western markets and its consumers. Led by a global vision and a drive to maintain quality upto international standards, our company strives to produce products conforming to international standards.

We would like more distributors to ADD their companies to our network and build beautiful Smiles all over!!

Cold Cure / Rapid Repair Powder & Liquid CC





Packing	Product Code
Cold Cure Powder Clear	(100 gms-6201/100)-(400 gms-6201/400)
Cold Cure Powder Pink	(100 gms-6202/100)-(400 gms-6202/400)
Cold Cure Powder Light Veined	(100 gms-6203/100)-(400 gms-6203/400)
Cold Cure Liquid	(100 ml-6301/100)-(400 ml-6301/400)
Cold Cure Universal Pack Clear	400 gms + 400 ml(6504/001)
Cold Cure Universal Pack Pink	400 gms + 400 ml(6505/001)
Cold Cure Universal Pack Light Veined	400 gms + 400 ml(6506/001)

Properties & Advantages:

- This is a 2 part rigid dental acrylic used for making dental appliances.
- >> Contains both powder and liquid components.
- ▶ Pink Dental Acrylic stimulates the fleshy gum tissue of the mouth.
- **▶** By combining this with pre-made acrylic teeth, or with your own castings made from tooth shade acrylic can make or repair dentures.
- For best results use with MAARC cold cure liquid.

Indication for use: Used for full or partial upper or lower dentures.

Shelf life and storage: Best before 36 months under proper storage.



Description & Features:

- ▶ Working time of 4 minutes
- ▶ Rapid bench cure time of 14 minutes
- ➤ Also useful in orthodontic work

Cold Mould Seal / Separating Medium (€ 6401/110

Properties & Advantages:

Packing

110 ml

5 Liter

- ▶ Cold Mould Seal is a specially developed formulation of separating solution for universal use, suitable for both packing / pressing technique & cast moulding.
- Makes tough film with elasticity which is unbreakable in pressure.

6401/005

→ Film does not chip off from plaster at pressure. High quality of the raw material guarantees optimum separation with both techniques.

Shelf life and storage: Best before 36 months under proper storage.









Packing	Product Code
Heat Cure Powder Clear	(100 gms-6001/100)-(450 gms-6001/450)
Heat Cure Powder Pink	(100 gms-6002/100)-(450 gms-6002/450)
Heat Cure Powder Light Veined	(100 gms-6003/100)-(450 gms-6003/450)
Heat Cure Liquid	(100 ml-6101/100)-(225 ml-6101/225)
Heat Cure Universal Pack Clear	450 gms + 225 ml + 110 ml CMSL(6501/001)
Heat Cure Universal Pack Pink	450 gms + 225 ml + 110 ml CMSL(6502/001)
Heat Cure Universal Pack Light Veined	450 gms + 225 ml + 110 ml CMSL(6503/001)

Properties & Advantages:

- >> This is a 2 part rigid dental acrylic used for making dental appliances.
- >> Contains both powder and liquid components.
- ➤ Light Veined Pink Dental Acrylic stimulates the fleshy gum tissue of the mouth.

 Indication for use: Used for full or partial upper or lower dentures

Description & Features

- >> Heat cured acrylic denture base material.
- >> Simple mixing, easy to process.
- >> Cross linking ensures superior quality

Shelf life and storage : Best before 36 months under proper storage. For Best results, Use with MAARC Heat Cure Liquid.



R2 Alloy (Ni-Cr Ceramic Alloy for PFM Bonding)





It comes in the form of rectangular pcs

Composition: Ni: 62%, Cr: 25%, Mo: 11%, Si 1%, Rest 1%.(Beryllium free)

Technical Data: Coeffecient of expansion (0-500°C): range from

13.9-14.1µm/m°C

Packing :	Product Code
500gm	7002 / 500

Properties & Advantages:

- >> Optimum castability and bonding strength in PFM restorations.
- >> It produces a uniform oxide film at porcelain metal interface ensuring strong chemical bonding.
- → High mechanical strength & stability on account of refined grain structure.
- ▶ Good burnishability.
- **▶** Brighter shine and finish.
- > Suitable for open flame casting.

Ni – Cr Alloy / White Metal for Cast Restorations & Crown & Bridge (Non Precious Alloy)

It is in the form of rectangular pellets

Nominal Composition: Ni - 68%, Cr - 15%, Cu - 15% and Rest - 2%

Packing :	Product Code
500gm	7001 / 500

Properties & Advantages:

- >> It is beryllium free alloy with excellent castability.
- >> It is reusable with addition of appx 50% fresh metal during recasting.
- >> It can be cast in open flame.
- >> It has got good burnish ability and is easy to finish.
- ▶ It can be used for making cast crown bridge and pontics and frames for hard resin crowns.









Color: Transparent

Packing :	Product Code
100ml glass bottle	8103 / 100

Properties & Advantages:

- ▶ It produces a homogenous coating on the die stone and plaster attributed to its viscosity.
- >> It seals, hardens and water proofs stone and plaster.
- ▶ It protects models against abrasion specifically at marginal area which is most prone to breakage.
- >> The film does not disassociate from the die while building wax pattern due to heating.
- **▶** It is rapid setting and thereby reduces working time.

Shelf life and storage: Best before 24 months under proper storage.



Color: Silver / Golden.

Packing :	Product Code
100ml glass bottle with 10ml applicator bottle+50ml Thinner	(Golden 8101 / 100)-(Silver 8102 / 100)

Properties & Advantages

- ▶It produces a homogenous coating of 5¼ in single application on the die for crown and bridge restoration thereby assuring accurate seating of restorations.
- ▶ The film does not disassociate from the die while building wax pattern due to heating.
- ▶ It is rapid setting and thereby reduces working time.

Shelf life and storage: Best before 24 months under proper storage.



Die Separator / Die Lubricant 🤇





It is the final product that would be applied to a die prior to the formation of the wax pattern. It is an ideal release agent for wax, porcelain, and composite materials.

Packing :	Product Code
100m l	3003 / 100

Shelf life and storage : Best before 24 months under proper storage.

Thickness: 0.25 mm **Type:** Golden and Silver

Packing	Product Code
Golden 10 Pcs	5502 / 010
Silver 10 Pcs	5501 / 010

Properties & Advantages:

- ➤ Re-inforcement Mesh is made from 316 grade stainless steel plated with gold or plain polished .
- >> These Mesh are very much flexible and can easily adapt any shape.
- >> It is 0.25 mm thin so can be easily layered between acrylic.
- **▶** It gives denture desirable strength.

















Thickness: 1.3mm - 1.5mm. (Custom Size / Thickness Catered),

Size: 160mm X 90 mm

Packing	Product Code
Modelling Wax Regular	(12 Sheet-2001/012)-(500 gms-2001/500)
Modelling Wax Red	(12 Sheet-2002/012)-(500 gms-2002/500)
Plasti wax Type 6	(12 Sheet-2015/012)
HARD WAX	(12 Sheet-2010/012)-(500 gms-2010/500)

Color: It comes in various shades such as pink, light pink, red.

Type: Soft, Medium, Hard & Plasti Wax.

Properties & Advantages:

- >> The wax is formable at slight increased temperature.
- Sufficient hard and resistant on cooling to prevent occlusion changes during trying with patient before making Final Denture to prevent failure of final denture.
- ▶ Easy & Clean carve on teeth.
- **▶** Bends without cracking.

Shelf life and storage:

Best before 60 months under proper storage.

Thickness: Upper / Lower Base Plate: 1.3mm, 1.7mm. 2.0mm, 2.5mm and 3.0mm (+/- 0.15mm) (Custom Thickness Catered)

Colour: Brown / Pink / Aluminium

Packing	Product Code
Brown Base Plate (8 Upper + 4 Lower)	1001/084
Yellow base Plate (8 Upper + 4 Lower)	1024/084

Type: It comes in different formulations for different climatic conditions (LMP / HMP / Al. Filled)

Properties & Advantage

- ▶ They are ready to use Maxillary & Mandibular plates 8 Upper & 4 Lower.
- **▶** Optimum heat resistance to minimize overheating of the Plate.
- **▶** Easily adaptable on heating and easy to work with.
- >> Retain the shape stability and do not wrap back on cooling.
- >> Optimum mechanical strength to resist breaking during Try-in denture
- Does not adhere to plaster.

Shelf life and storage: Best before 24 months. Store below 30°C under proper storage and avoid from direct sunlight.

Shellac Base Plates (



Dental Stone Plaster (Green, Yellow & Off White)





For fabricating custom Trays, mouth guards and vaccum formed splints. For denture repairs. For studycasts. For lab procedures such as articulator cast mounting. For processing full or partial acrylic dentures.

Packing :	Product Code
(Yellow)	
30kg polyester gunny / paper bag.	4001 / 030
(Green)	
30kg polyester gunny / paper bag.	4002 / 030
(Off White)	
30kg polyester gunny / paper bag.	4005 / 030

Type: III - Specification:

Water/Powder	30m l /100gm
MixingTime	60 seconds wet
WorkingTime	6 minutes dry
SettingTime	8 to 12 minutes

Hard Set (Initial)	30 Minutes
Wet Compressive Strength	250 kg/cm2
Dry Compressive Strength	500 kg/cm2
Setting Expansion (Max)	0.08%

Shelf life and storage:

Best before 24 months when stored in cool and dry place.

For crown and bridge, dies and cast partials. For basing dental arches and opposing models. For making study models. In use where economical yet tough, non brittle and durable models are required.

Packing :	Product Code
(White) 30kg polyester gunny / paper bag.	4102 / 030

Type: IV - Specification:

Water/ Powder	24-26m l /100gm
MixingTime	45 seconds wet
WorkingTime	5 minutes dry
SettingTime	8-12 minutes

Hard Set (Initial)	30 Minutes
Wet Compressive Strength	400 kg/cm2
Dry Compressive Strength	800 kg/cm2
Linear Expansion (Max)	0.08%

Shelf life and storage: Best before 24 months when stored in cool and dry place.

Ortho White Stone Plaster CC









For crown and bridge, dies and cast partials. For basing dental arches and opposing models. For making study models. In use where economical yet tough, non brittle and durable models are required.

Packing :	Product Code
30kg polyester gunny / paper bag.	4003 / 030

Type: III - Specification:

D98.4x22mm

D98.4x25mm

<i>7</i> 1	
Water/ Powder	26-28ml/100gm
MixingTime	60 seconds wet
WorkingTime	5 minutes dry
SettingTime	8-10 minutes


Hard Set (Initial)	30 Minutes
Setting Expansion	0.2%
Compressive 1h.Mpa	35
Compressive 24h.Mpa	75

Shelf life and storage:

Best before 24 months when stored in cool and dry place.

Girrbach System Zirkonzahn System **Open System** Other D98x10mm D-wax 13mm D95x10mm D98.4x10mm D98x12mm D98.4x12mm D-wax 20mm D95x12mm D98x14mm D95x14mm D98.4x14mm D98x16mm D95x16mm D98.4x16mm D98x18mm D95x18mm D98.4x18mm D98x20mm D95x20mm D98.4x20mm

Milling Wax Blanks (

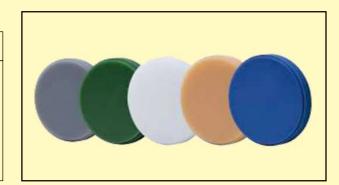

Packing: 1 Pc Per Pack

D98x22mm

D98x25mm

Colours Available: White, Ivory, Blue, Pink, Green, Tooth Colour, Grey etc.

D95x22mm

D95x25mm

Sprue / Wire wax (€





Packing :	Product Code
1.5mm	2301 / 500
2.0mm	2302 / 500
2.5mm	2303 / 500
3.0mm	2304 / 500
4.0mm	2305 / 500
5.0mm	2306 / 500

Thickness: It is in the wire form available in diameter of 1.5mm, 2mm, 2.5mm, 3mm, 4mm, 5mm and 6mm (Custom Thickness / Shapes Catered)

Packing: 500 gms (Can be customized)

Properties & Advantages:

- >> Stiff enough to hold the wax patterns and does not deform in shape.
- ➤ Easy burn out without residue

Shelf life and storage: Best before 60 months under proper storage.

Description: Pre-formed lingual connector bars. Ideal for clasps, flange and lingual connector areas. Used by dental labs while waxing a denture.

Thickness: These are wax strips - 160mm long x 15mm high, with taper across 15mm - one edge - 2.25mm and another 1.1mm.

Packing :	Product Code
453gms (1lb)	2205 / 453

Shelf life and storage: Best before 60 months under proper storage.

Lingual Bar Wax CE









Description: Utility Wax are in form of ropes / wires, easily shaped waxes that are specially suited for fill-ins and wax pattern repair. It can also be used to easily correct mistakes by "blending-in" with your original wax part. Extremely pliable at room temperature. It can also be used to fill sprue bases.

Thickness: 5.0 diameter round x 30cm long

Packing :	Product Code
12 Strips	2308 / 012
250 gms	2308 / 250

Shelf Life: Best before 60 months under proper storage.




Which retention profile shape of the various MAARC wax retention is determined by the specific indication and clinical requirement. Excellent restoration result can be achieved in both mandibular and maxillary modelling. Both partial and total dental prosthesis can be modelled easily and efficiently.

MAARC wax retentions are available in packs of 10 pieces.

Shelf Life: Best before 24 months under proper storage.



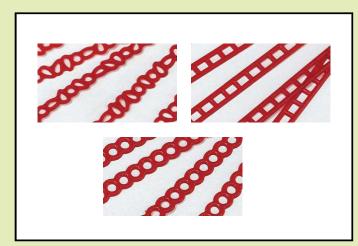







Packing :	Product Code
10 pcs	2952/010

Packing: 10 Pcs

Slender shape with high stability is one of the requirements met superbly by these retentions. These retentions are characterised by two

Important Aspects: the complexity of wax-up is reduced, thus saving time, and the partial denture prosthesis acquires optimal stability and wearing comfort. The grid retentions with holes can also be integrated very well as reinforcing grids in total maxillary dentures made of plastic.

Shelf Life: Best before 60 months under proper storage.

Packing	Product Code			
	2901/010	2904/010	2907/010	
10 pcs	2902/010	2905/010	2908/010	2910/010
	2903/010	2906/010	2909/010	

Packing: 10 Pcs

Wax clasp profiles for molars and premolars are designed to ensure good hold and stability for dentures. Half-teardop shaped cross-sections are effective at preventing food residues from getting stuck. All clasp profiles have an optimal anatomical shape and are made of easily adapted wax. Adhesion to the investment material model is optimal.

MAARC wax clasp profiles help save time during modelling, and highly customised adaptation of the wax shapes can be effected by shortening or lengthening them.

Shelf Life: Best before 60 months under proper storage.

Clasp wax (€





Inlay / Onlay Pattern Wax (



Types :	Product Code	Description
Cervical Wax (1) Orange	2401 / 200	Soft type with degree of fluidity that will allow it to flow and fit perfectly in cervical areas, high adaptation and stability to surfaces. It will wet model surface without shrinkage.
Cervical Wax Wax (2) Red/Blue	2402 / 200	Harder than orange wax.
Cervical Wax (3) Blue	2403 / 200	Standard inlay wax wherein the pattern comes out easily from the plaster surface.
Undercut Wax (4) Green	2404 / 200	Harder than red wax used where hard pattern waxes are preferred.
Undercut Wax (5) Violet/Ivory	2405 / 200	Harder than green wax - used for coping of patterns to add mechanical strength and as undercut wax.
Sculpturing Wax (6) Coffee	2406 / 200	Used for undercuts and adding mechanical stability to wax patterns of jacket crowns.
Carving Wax (7) Magenta	2407 / 200	Hard inlay wax. With medium flow for built up.
Carvical Wax (8) Ivory Inlay wax	2408 / 200	A white, enamel-like wax used to makes a more aesthetic case presentation (wax-up). Easy carving properties. Also aids in the fabrication of temporaries.

Packing: 200 gms

Properties & Advantages:

- ➤ The basic classification from our side is being done based on hardness and stiffness & flow of waxes.
- ➤ Softer waxes with high flow are more flexible and are often used.

Shelf Life: Best before 60 months under proper storage.

Die Stone Plaster (Pink) 🥨





For crown and bridge, dies and cast partials. For basing dental arches and opposing models. For making study models. In use where economical yet tough, non brittle and durable models are required.

Packing	Product Code
30 kg	4101/030

Type: IV - Specification:

Water/ Powder	24-26ml/100gm
MixingTime	45 seconds wet
WorkingTime	5 minutes dry
SettingTime	8-12 minutes

Hard Set (Initial)	30 Minutes
Wet Compressive Strength	400 kg/cm2
Dry Compressive Strength	800 kg/cm2
Linear Expansion (Max)	0.08%

Shelf life and storage:

Best before 24 months when stored in cool and dry place.

Packing Product Code Packing 250 gms 2206/250 Packing 200 gms 2206/200

Thickness: 0.3 to 0.6 mm Size: 160mm x 90 mm

Colour: 'Maarc' Casting Wax sheets are thin wax sheets. Available in

Green color for cobalt chrome cast partial dentures.

Packing: Available in 200gms pack. The Sheets are wrapped in non adhesive Silicone Paper to Facilitate easy removal.

Properties & Advantages:

- ► The Sheets are sufficiently thin to have thinner casting as possible without defects to reduce machining.
- >> They easily burn out during castings.

Shelf life and storage:

Best before 60 months under proper storage.

Casting Wax Sheets @



INSTITUTIONAL RANGE







Properties & Advantages:

- → The Wax Blocks are easy to carve.
- ➤ The special feature being the wax is sufficiently hard that it does not chip off at edges on applying pressure.
- ▶ Free from air bubbles.

Shelf life and storage:

Best before 60 months under proper storage.

Packing	Prod	uct Code
40 Blocks	White	2501/040
	Red	2502/040
	Blue	2503/040
	Yellow	2504/040

'MAARC' Carving wax comes in small blocks slugs and in various colors like Translucent White, Yellow, Red, Blue

Thickness: 38 x 13 x 13 mm

Packing Product Code 12 strips 2309/012

Thickness: 6mm x 6mm x 250mm appx.

Packing: It comes in the pack of 12 strips wrapped in non adhesive Silicone Paper to facilitate easy removal.

Properties & Advantages:

- ▶ It is sticky enough to easily adapt at the periphery or edge of alginate impression / final impression by hand pressure.
- ▶ It creates the base of the cast and preserves the outer portion of the cast.

Shelf life and storage: Best before 60 months under proper storage.





INSTITUTIONAL RANGE

Green Tracing Sticks (





Thickness: 113 x 9 mm dia weighing appx. 10gms.

	J J 1 1
Packing	Product Code
10 Sticks	1501/010
15 Sticks	1501/015
100 Sticks	1501/100

Indication:

- ➤ Low fusing material primarily used for border extensions on impression trays.
- → High viscosity also makes it ideal when recording the full depth of the sulcus.
- ➤ Modify dentate tray by building up material where tooth has been extracted.
- Make functional and relining impressions.

Properties & Advantages:

- >> Suitable for muscle trimming work.
- >> Fast and uniform plasticity.
- ▶ Quick hardening in cold water.
- Wide working area.

Shelf life and storage: Best before 60 months under proper storage.





Thickness: L 114 x 8 mm dia of Big Sticks

Packing	Product Code
10 Sticks (70 gms)	2403/010

Properties & Advantages:

- >> It comes in stick form.
- **▶** It is used in institutes by students for easy application.

Shelf life and storage:

Best before 60 months under proper storage.





Thickness: 114 x 8 mm dia sticks weighing appx. 6gms. / Pc

Packing	Product Code
10 pcs	2701/010
200 gms	2701/200


Properties & Advantages:

- ▶ It is blending of various waxes and resins and it is yellow in color.
- >> It melts when heated and adheres to the surface over which it is applied.
- >> It holds the components prior to fixation and gives a clean burnout.

Shelf life and storage: Best before 60 months under proper storage.

INSTITUTIONAL RANGE

Spacer Wax (6

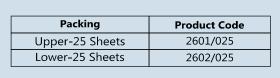



Thickness: 0.4mm (+/-0.1mm)

Packing:

Upper and Lower Spacer Wax are packed separately. Wrapped in non adhesive Silicone Paper to facilitate easy removal.

Properties & Advantages:

- → Thin wax sheets of 0.3mm thickness shape of Upper and Lower Base Plates.
- >> It increases the thickness of wax rims while making dentures using base plates.
- >> Provide uniform thickness of impression material.
- >> Use of spacer wax prevents direct contact of acrylic trays with soft tissues.

Shelf life and storage:

Best before 60 months under proper storage.

Packing	Product Code
6 inch-6 Sheets	2203/006
12 inch-6 Sheet	s 2204/006

Thickness: 160 x 45 x 1.5mm

Packing : Pack of 6 Sheets wrapped in silicon paper to facilitate easy removal.

Properties & Advantages:

- **▶** It is soft and tacky for boxing in impressions.
- ▶ Adheres easily to most of the surfaces.

Shelf life and storage:

Best before 60 months under proper storage.





Proudly Exporting to 40+ Countries Worldwide





Shiva Products

ISO 13485: 2016 (TUV SUD Certified)

Manufactured by: Shiva Products

H.O.: 214, Bluebell Indl. Estate, Chinchpada, Vasai (East), Dist. Palghar 401 208, Maharashtra, India.

∴ B-18, Dwarka Indl.Est., Naik Pada, Waliv, Vasai (East), Dist. Palghar 401 208, Maharashtra, India.

Branch Office Shop No. 7, Gola Market, Pratap Street, Nr. Golcha Cinema, Darya Ganj, New Dehli-110 002. (INDIA)





