

Global Market



i3 Automated External Defibrillator

Compliance with ERC/AHA guidelines

Amoul® Ambulanc (Shenzhen) Tech. Co., Ltd.



Tel: +86 755-26072286 E-mail: info@ambu-lanc.com
Fax: +86 755-23016012 website: www.ambulgroup.com

Address: Building 5C, Skyworth Innovation Industry Park, Tangtou 1st Road, Baoan District, Shenzhen, P.R.China

© Copyright 2021 Ambulanc(Shenzhen) Technology Co., Ltd
SH SZAB I3 AED/CE20210927-V1.0 Printed in China

Amoul®
We Innovate to Save Lives

Automated External Defibrillator

can save life when a heart arrest happens

i3



Real-time battery status



Voice and visual instructions guide the rescuer through the entire defibrillation process



Smart Biphasic Truncated Exponential shock waveform



Remotely management



Portable size & Lightweight (2.3kg)




1 Turn on AED



2 Place the pads on the patient as directed.



3 Follow the voice instructions and press the shock button

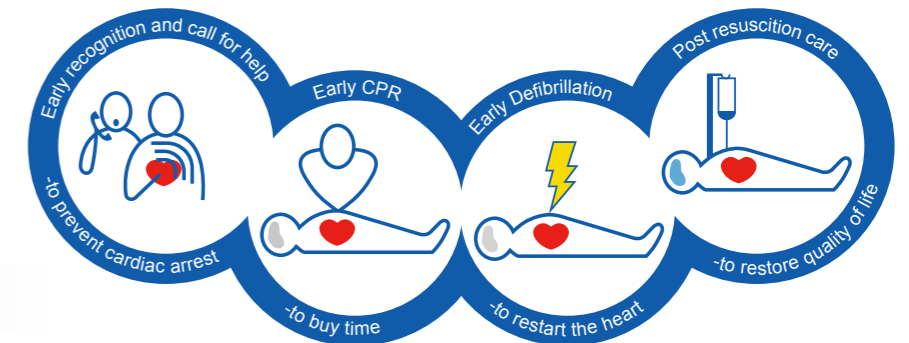

Features

- Biphasic Truncated Exponential (BTE) Shock Waveform
- High capacity and long-life lithium manganese dioxide battery, minimum 200 shocks, 5 years valid
- Visualized status of battery and system
- Detailed voice instructions and illustrations makes rescue easier
- Automated tracking and direct data transfer via built-in 4G module (optional)
- Data storage of 24hrs ECG and up to 1000 events
- Deliver a 200J shock to adults and 50J shock to pediatric
- Daily self-tests, report device status in time, ideal for remote management.